



# RetractOrtho

## Retract, Reduce, Repair

---



Metaphyseal System shown above

**The RetractOrtho Retractor System is a radiolucent, self-retaining surgical retractor that enables efficient reduction and repair of metaphyseal and diaphyseal bone fractures. The system is delivered sterile and ready for single patient use.**

### Hands-free

Once set-up, hands-free retraction enables the surgeon to work independently, optimizing workflow efficiency and visualization.

### Self-retaining




Self-retaining design eliminates the need for ongoing retractor management and reduces the need to repeatedly insert, adjust and remove retractor, minimizing soft tissue injury.

### Stable

Textured backing plate provides a stable platform for sharp forceps to aid fracture reduction.

### Radiolucent

Radiolucent system allows for visualization of surgical field under fluoroscopic view without interference of metal retractors. Implant positioning is more predictable, improving surgical efficiency and workflow.

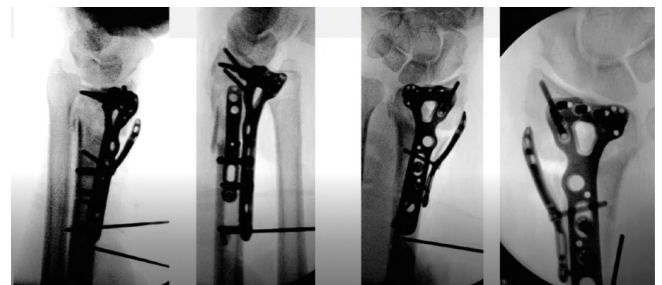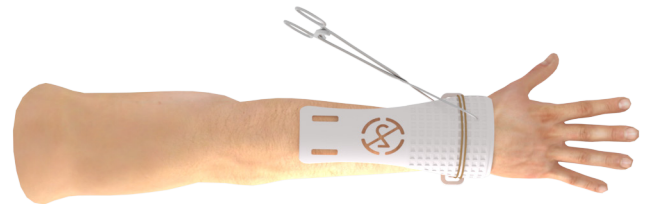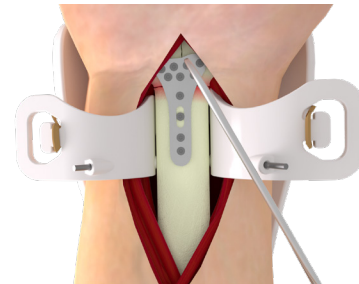


**To order, contact Vital Hospital Systems:**


**Phone: (314) 394-0242 • Fax: (314) 394-1543 • Email: [customercare@vitalhospital.com](mailto:customercare@vitalhospital.com)**

# RetractOrtho Retractor System Components

## Metaphyseal Retractor System Components

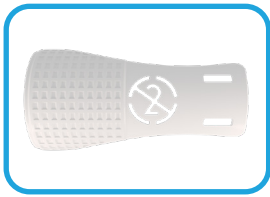


Textured Backing Plate

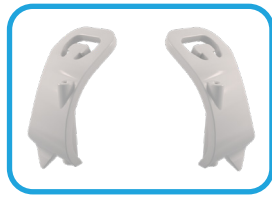


Asymmetric Single Pin Retractors (1 ea)



Universal Pin and Place Guide

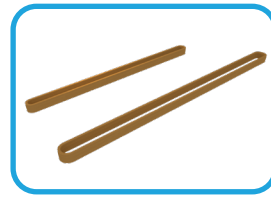


4" and 6" Elastic Bands (1 ea)

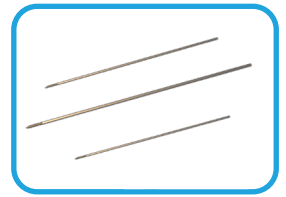


Three (3) Custom Threaded K-Wires


## Diaphyseal Retractor System Components

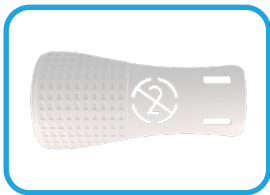


Textured Backing Plate

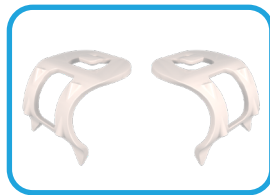

Symmetric Double Pin Retractors (2)

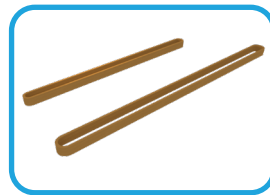

4" and 6" Elastic Bands (1 ea)

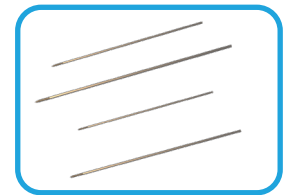

Four (4) Custom Threaded K-Wires

Note: Surgical Reduction Clamps are not included

## Ordering Information

Manufacturer Product Code	Quantity	Description
R3-META-M-01	1	Retractor kit for metaphyseal surgical procedures. Single unit.
R3-DIA-M-01	1	Retractor kit for diaphyseal surgical procedures. Single unit.
R3-META-M-05	5	Retractor kit for metaphyseal surgical procedures. Five pack.
R3-DIA-M-05	5	Retractor kit for diaphyseal surgical procedures. Five pack.
R3-PIN-U-01	1	Universal Pin & Place Guide. Single unit.
R3-PIN-U-05	5	Universal Pin & Place Guide. Five pack.